

SPEAKING


It is easy to control your speech with only one hand.

NO TOOLS REQUIRED.

Adjustable top head: Twist for fine-tuning of the sound, if needed. **WARNING** - over adjusting of the head will result in diminished tone quality; be careful.

On/Off Button: To speak, push the button and release when you are finished.

Pitch Wheel: (on the left) Turn the wheel for easy adjustment to a higher or lower pitch. No need to further adjust after you have achieved your desired pitch level.

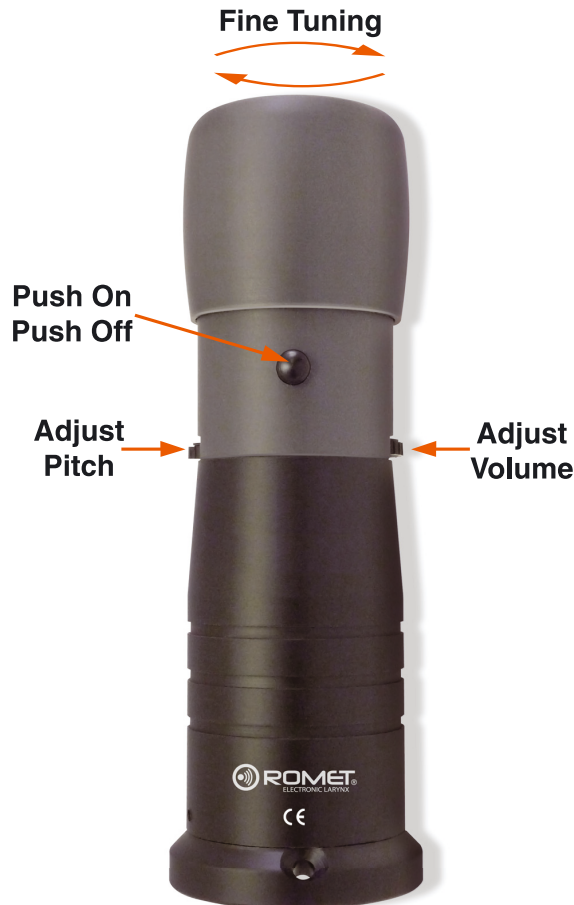

Volume Wheel: (on the right) **Note:** unit is shipped with the volume **OFF**. The wheel must first be turned up to activate volume; then used for easy adjustment while speaking to increase or decrease your speaking volume for emphasis.

Practice: Place the Romet® R355 Electronic Larynx to your neck and practice finding your "sweet spot" to achieve your best vocal tone.



NOTE: This spot can be many places on either your neck or the side of your mouth; each person is different depending upon your skin texture and the amount of radiation you have received. Look into a mirror and repeat simple sentences to achieve vocal clarity. If possible, record your voice and play it back to hear how you sound. Practice regularly!

CHARGING

Charging your unit couldn't be easier: Simply plug in your charging pad and place the R355 unit **centered** on the pad (as shown). While charging is in progress, the light will flash red. When charging is complete, the light will turn green. If the light does not turn on when you place the unit on the pad, it still has a sufficient charge. The unit should be fully charged within approximately 3 hours.



ORAL ADAPTER INSTRUCTIONS

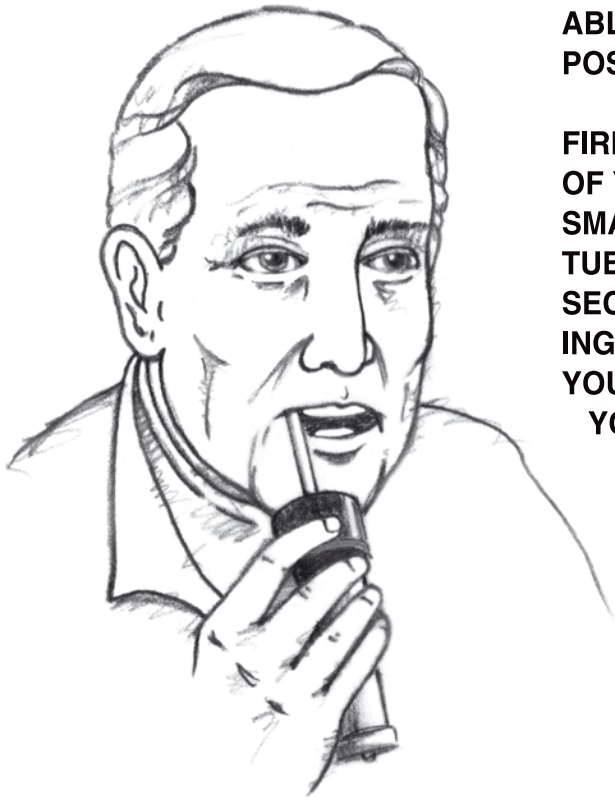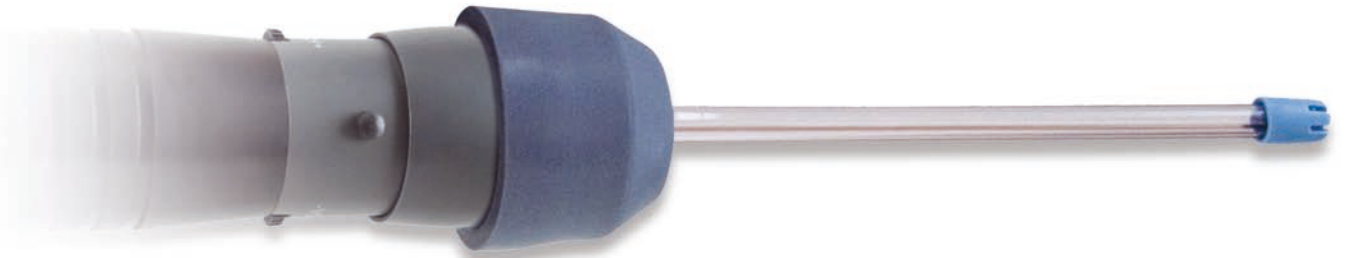

THE ORAL ADAPTER AND TUBES PROVIDED WITH YOUR ROMET ELECTRONIC LARYNX WILL ENABLE YOU TO BE ABLE TO BEGIN SPEAKING WHILE YOUR NECK IS HEALING POST SURGERY.

FIRMLY PLACE THE RUBBER CAP ADAPTER ON THE TOP OF YOUR ELECTRONIC LARYNX. INSERT A TUBE INTO THE SMALL HOLE ON THE TOP OF THE ADAPTER. PUSH THE TUBE FIRMLY INTO THE HOLE. THE TUBE MUST BE SECURE IN THE ADAPTER HOLE WHILE YOU ARE SPEAKING. PLACE THE END OF THE TUBE INTO THE CORNER OF YOUR MOUTH POINTING THE TIP TOWARDS THE TOP OF YOUR MOUTH.

ACTIVATE THE ON-OFF BUTTON. PRACTICE YOUR SPEECH WITH SIMPLE VOWELS OR SHORT SOUNDS. TRY A FEW SIMPLE SENTENCES.

IF POSSIBLE, SPEAK SLOWLY INTO A TAPE RECORDER. REPLAY YOUR SPEECH.

CLEAR TUBES ARE FITTED WITH TIPS TO PREVENT ANY FOOD PARTICLES FROM LODGING INTO THE TUBE WHILE BEING USED. YOU CAN EASILY REMOVE THIS TIP OR SHORTEN THE LENGTH OF THE TUBE WITH SCISSORS.

AFTER FINISHED SPEAKING, RINSE THE TUBES WITH COOL CLEAR WATER. BE SURE TO CLEAN THE TOP OF YOUR ELECTRONIC LARYNX TO REMOVE ANY SALIVA RESIDUE OR FOOD PARTICIPLES.



SANITAS GRUP Ltd.

1504 Sofia, Bulgaria, 25 "Tzar Osvoboditel" bul., ap. 1
office number: +359 2 943-38-38, office/fax: +359 2 944-16-50